


Surgical Aspirators – fits 11mm High Volume Valve



overall length 6.25 inches

Surgical aspirator made of 1/4" diameter stainless steel. Tapered tip prevents clogging. Larger tube size allows for more volume of suction. Recommended cleaning brush size (Pg. 25): Medium.

Opening	Fits Into	Product #	Price
1.5 mm	HVE Cut-Off Valve	46P1A	\$30.00
2.5 mm	HVE Cut-Off Valve	46P2A	\$30.00
3.0 mm	HVE Cut-Off Valve	46P3A	\$30.00

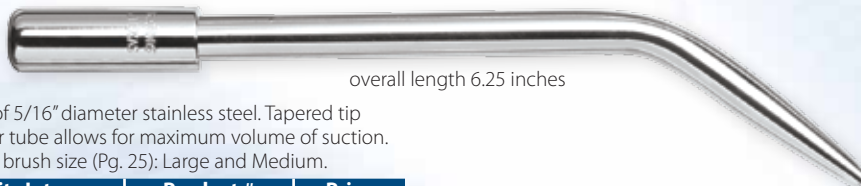


overall length 6.25 inches

Surgical aspirator made of 5/16" diameter stainless steel. Tapered tip prevents clogging. Larger tube allows for maximum volume of suction. Recommended cleaning brush size (Pg. 25): Large and Medium.

Opening	Fits Into	Product #	Price
2.5 mm	HVE Cut-Off Valve	56P2	\$30.00
3.0 mm	HVE Cut-Off Valve	56P3	\$30.00
4.0 mm	HVE Cut-Off Valve	56P4	\$30.00

- 1.5 mm ●
- 2.5 mm ●
- 3.0 mm ●
- 4.0 mm ●

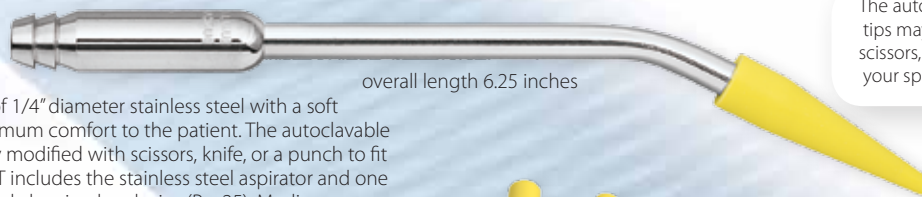


overall length 6.25 inches


Surgical aspirator, made of 1/4" diameter stainless steel with a soft silicone tip provides maximum comfort to the patient. The autoclavable silicone tips may be easily modified with scissors, knife, or a punch to fit your specific needs. 46YST includes the stainless steel aspirator and one silicone tip. Recommended cleaning brush size (Pg. 25): Medium.

Opening	Fits Into	Product #	Price
Variable	HVE Cut-Off Valve	46YST	\$30.00

The autoclavable silicone tips may be easily modified with scissors, knife, or a punch to fit your specific needs.

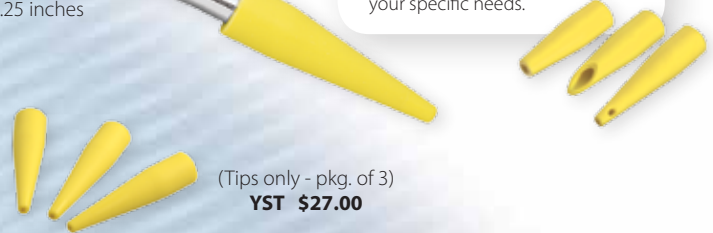

(Tips only - pkg. of 3)
YST \$27.00

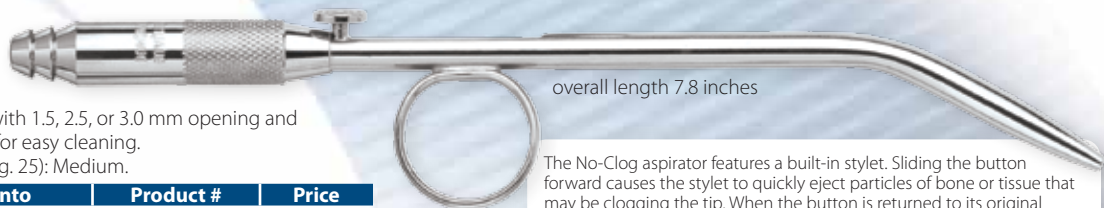

overall length 7.8 inches


No-Clog surgical aspirators with 1.5, 2.5, or 3.0 mm opening and built-in stylet. Disassembles for easy cleaning. Recommended brush size (Pg. 25): Medium.

Opening	Fits Into	Product #	Price
1.5 mm	HVE Cut-Off Valve	B1	\$105.00
2.5 mm	HVE Cut-Off Valve	B2	\$105.00
3.0 mm	HVE Cut-Off Valve	B3	\$105.00

The No-Clog aspirator features a built-in stylet. Sliding the button forward causes the stylet to quickly eject particles of bone or tissue that may be clogging the tip. When the button is returned to its original position, suction resumes. The adapter end will fit high volume cut-off valves or surgical tubing.



Titanium Surgical Aspirators – fits 11mm High Volume Valve

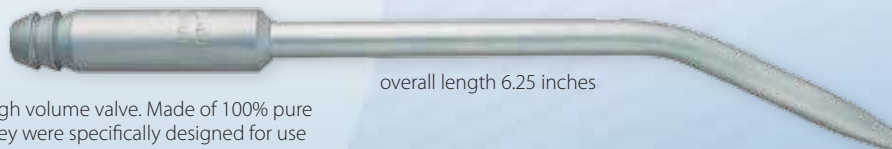

overall length 6.25 inches

Titanium aspirators for use with high volume valve. Made of 100% pure Titanium, not coated or plated. They were specifically designed for use with Titanium implants. Recommended brush size (Pg. 25): Medium.

Opening	Fits Into	Product #	Price
2.5 mm	HVE Cut-Off Valve	T-46P2A	\$125.00
3.0 mm	HVE Cut-Off Valve	T-46P3A	\$125.00